

## Certification



Promoting member of ISMST



CE

## Specification

Interface	10.4" LCD touch screen
Pressure Range	1.0 – 4.0 bar
Frequency Range	1 – 17 Hz
Pulse	1-15000
Dimension	52.5(L) × 47.6(W) × 110(H) cm
Standard Weight	27.5 kg
Preset Protocols	History Download / Yes
User-defined Protocols	Yes


# Longest

Guangzhou Longest Science & Technology Co., Ltd.

Address: #96 Chuangqiang Road, Ningxi Street, Zengcheng District, Guangzhou, China 511399

Tel: 86-20-66353999 Fax: 86-20-66353920

E-mail: [export@longest.cn](mailto:export@longest.cn) Website: [www.gzlongest.com](http://www.gzlongest.com)

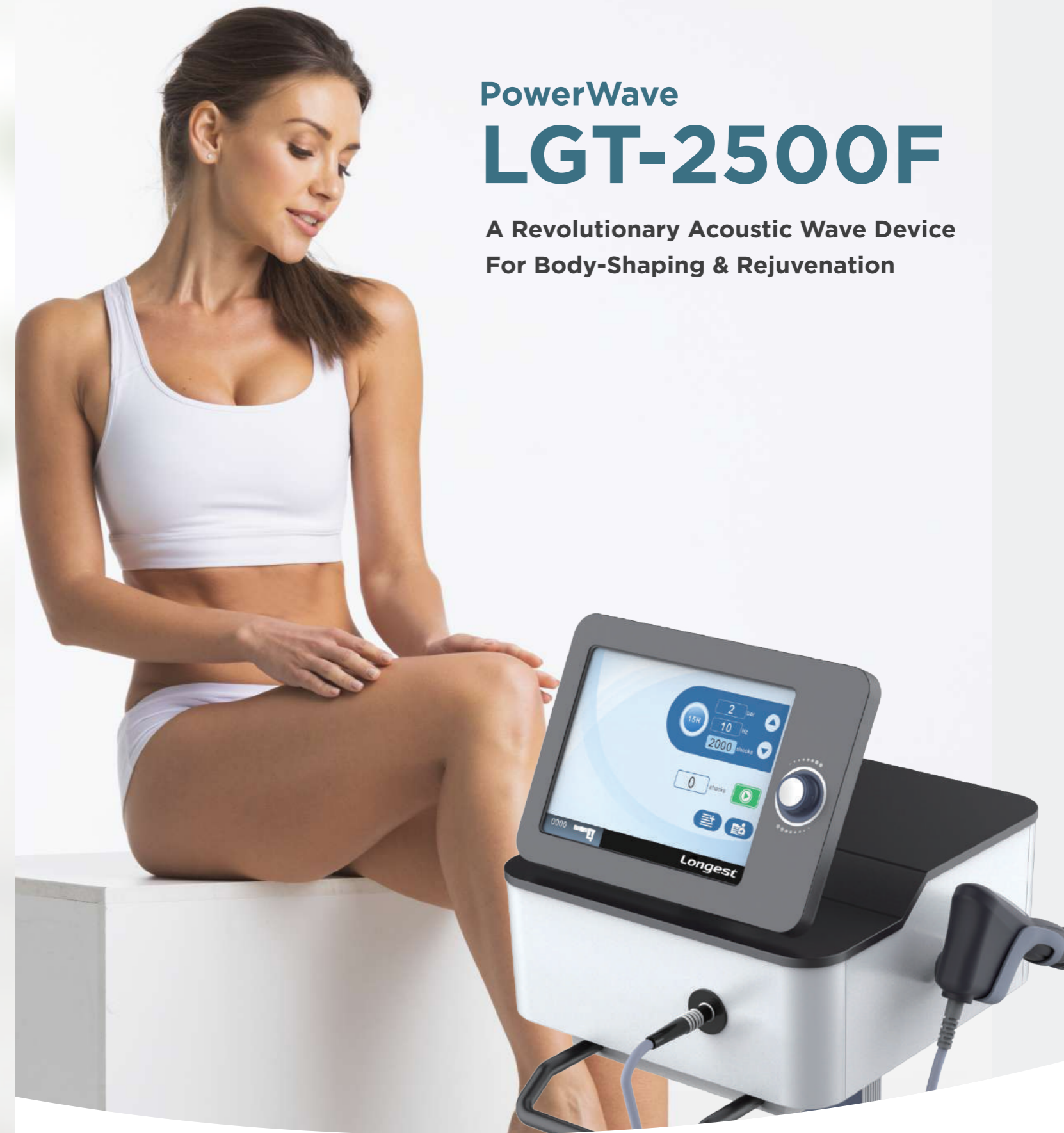

## PowerWave LGT-2500F

A Revolutionary Acoustic Wave Device  
For Body-Shaping & Rejuvenation

Health Life Longest Care

# Longest



## LGT-2500F High-performance Acoustic Wave Therapy Machine

LGT-2500F is a professional radial shockwave therapy machine, which can effectively improve skin texture, enhance elasticity, shape body, firm tissue, smooth scars, and eliminate cellulite.

This acoustic wave therapy machine combines acoustic sound waves with the body's own ability to repair tissue by stimulating circulation in all layers of the skin, blood vessels, connective tissue, and muscles.



### Features

- ⊗ **Foldable Touch Screen**  
10.4" LCD color screen; clear and intuitive display
- ⊗ **Rotary Knob**  
Quick setting and sensitive response facilitate operation
- ⊗ **Integrated Air Compressor**  
German KNF technology; silent and stable output
- ⊗ **Built-in Air Tank**  
Ensuring the stable, comfortable and more powerful output
- ⊗ **Handlebar**  
Holding design facilitates movement of the device
- ⊗ **Ergonomic Handpiece**  
Ergonomic grip and elimination of backward; 2 million shocks tested projectile
- ⊗ **Metal Shell**  
Fashionable design and high quality
- ⊗ **360° Lockable Wheel**  
Silent and flexible movement; lockable design anchors the device when it's placed

### Application



#### FACE

Firming  
Anti-ageing  
Redensification  
Wrinkles removal  
Tonification



#### BODY

Body-shaping  
Anti-cellulite  
Localised fat  
Flaccidity  
Fluid retention  
Stretch marks



#### MEN'S CARE

Erectile dysfunction  
Peyronie's disease  
Chronic prostatitis  
Benign prostatic hyperplasia

## High Profitability - Lease Mode

LGT-2500F supports leasing, and this function accomplishes recirculation of the device. Owners can lend it to end users or clinics to make more profit. Each lease period is 30 days. Costumers can renew a period by inputting the PIN code provided by the supplier.



## Why Choose to Lease?

Acquiring medical devices requires a considerable investment to remain competitive. With a simple and fixed monthly payment plan, a lease can enable a business to match scheduled payments with a practice's anticipated revenue stream.

- ✓ Avoid delays and simplify the complex process associated with the standard loan application process
- ✓ Structured fixed payments over a short or long-term basis
- ✓ Reduce budget in a lab experimentl
- ✓ Leasing does not impact a business credit line
- ✓ Lease and conserve valuable capital

## Advantages

- ✓ **Results**  
Lasting results for patients seeking no-drug, non-surgical solutions.
- ✓ **High ROI**  
In-office treatment offers a more profitable solution without the long downtime of complex surgeries.
- ✓ **Simplicity**  
Quick, easy-to-perform, delegable procedure.
- ✓ **Clinical Safety**  
Non-invasive, no recovery time, drug-free and painless.



## How It Works

Pulse Wave Therapy utilizes an internal ballistic projectile to transfer pulse wave energy through the handpiece. An acoustic sound wave, which diffuses energy as it travels away from the applicator tip, creates a therapeutic effect as it passes through the target tissue.

## Profitable Procedure for Practice Owners

- ✓ In-office Clinical Support
- ✓ 1 Year Limited Warranty on Hand Pieces
- ✓ 2 Year Warranty on Console
- ✓ Easy Return on Investment
- ✓ Low Acquisition Cost
- ✓ Low Maintenance Cost
- ✓ High Patient Satisfaction

## Standard Accessories

- Cleaning Brush
- Handpiece (with 15mm Radial Transmitter)
- Metal Pipe
- Screwdriver
- Spare Accessory Kit
- Projectile
- 36mm Radial Transmitter



## Optional Transmitters

Transmitters						
<b>Type</b>	6 A ( D = 6 mm Acupuncture )*	15 R ( D = 15 mm Radial )	15 F ( D = 15 mm Focus )*	15 D ( D = 15 mm Deep )*	20 R ( D = 20 mm Radial )*	36 R ( D = 36 mm Radial )
<b>Penetration Depth</b>	0 – 40 mm	0 – 40 mm	0 – 20 mm	0 – 60 mm	0 – 50 mm	0 – 50 mm

Notice: " \* " Stands for Optional Accessories.